

Vacumed® Double wall Sodium Citrate test tube

FL Medical works to satisfy customers' needs, The company is continuously updating to improve the products and realize new projects.
At the basis of this item there are FL Medical values

Innovation
Quality
Efficiency

The main characteristics of the Vacumed® Double wall Sodium Citrate test tube are

INNER TUBE

Material: Polypropylene

Capacity: 1,8 ml or 2,7 ml

**The Polypropylene permits a better sodium citrate maintenance
guaranteeing the correct additive/blood volume ratio**

OUTER TUBE

Dimension: 13 x 75 mm

Material: PET

**PET retains vacuum ensuring the blood
volume/additive ratio**

CAP


Dimension: 12 mm

Shield Material: Polipropylene

Stopper Material: Rubber

Colour: Light Blue

**The main features of the cap are the perfect design and the high quality of
materials which provide a perfect tightness**

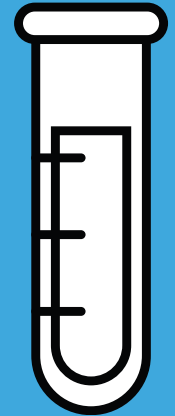

Code	Sodium Citrate	Aspiration Volume	Packaging
42450	3,2%	1,8 ml	100 pcs per rack 1.200 pcs per box
42451	3,2%	2,7 ml	100 pcs per rack 1.200 pcs per box
42490	3,8%	1,8 ml	100 pcs per rack 1.200 pcs per box
42491	3,8%	2,7 ml	100 pcs per rack 1.200 pcs per box